
seca 834

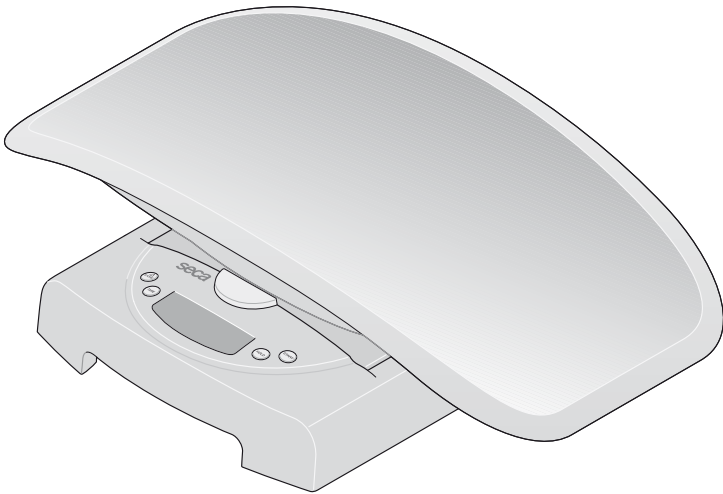D	<i>Bedienungsanleitung und Garantieerklärung</i>	3
GB	<i>Instruction manual and guarantee</i>	11
F	<i>Mode d'emploi et garantie</i>	19
I	<i>Manuale di istruzioni e garanzia</i>	27
E	<i>Manual de instrucciones y garantía</i>	35
DK	<i>Betjeningsvejledning og garantibevis</i>	43
S	<i>Bruksanvisning och garanti</i>	51
N	<i>Bruksanvisning og garantierklæring</i>	59
FIN	<i>Käyttöohje ja takuu</i>	67
NL	<i>Bedieningshandleiding en garantieverklaring</i>	75
P	<i>Instruções de utilização e declaração de garantia</i>	83
GR	<i>Οδηγίες χειρισμού και εγγύηση</i>	91

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Full certification

In seca products you are purchasing not only technology matured over a century, but also quality confirmed by authorities, the legal system and by institutes.

seca products comply with European directives, standards and national laws. With seca, you are buying the future.

M

The products in these operating instructions comply with German legislation governing medical products and therefore with Directive 93/42/EEC of the Council of the European Community which is embodied in national laws throughout Europe.

Scales carrying this symbol are approved for medical purposes in the European Community. seca scales with this symbol meet the stringent quality and technical requirements placed on scales which can be calibrated.

Scales carrying this symbol are calibrated in accordance with Precision Class IIII of EC Directive 90/384/EEC.

Products carrying this symbol comply with the following directives and standards:

1. Directive 90/384/EEC governing non-automatic weighing instruments
2. Directive 93/42/EEC governing medical devices
3. DIN EN 45501 governing metrological aspects of non-automatic scales

seca's professionalism has also been officially recognized. The certificate from TÜV Product Service, the body responsible for medical devices, confirms that as a manufacturer of medical devices, seca rigorously complies with stringent legal requirements. seca's quality assurance system covers the development, production, sales and service of medical scales and measuring systems and ergometry.

seca helps the environment.

We are concerned to save natural resources. We therefore strive to save packaging material where meaningful. What is left over can conveniently be disposed of locally via the German Dual System (recycling programme).

1. Congratulations!

In the **seca 834** electronic baby scale you have acquired a highly-accurate yet simultaneously sturdy piece of equipment.

seca has been putting its experience at the service of health for over 150 years now, and as market leader in many countries of the world, is always setting new standards with its innovative developments for weighing and measuring.

The **seca 834** electronic baby scale is designed for use in a doctor's surgery or hospital. It serves both to establish birth weight and to check weight gain during the growth process.

The scale has Class III approval.

2. Safety


- Follow the instructions in the instructions for use.
- Keep the operating instructions and the declaration of conformity in a safe place.
- Never leave a baby lying on the scale unsupervised.
- Ensure that the scale is standing firmly on a smooth, level surface.
- Do not drop the scale or subject it to violent shocks.
- Use only the type of battery stated (see "Power supply" on page 12).
- Have scale serviced and re-calibrated on a regular basis (see "Servicing/recalibration" on page 16).
- Have repairs carried out only by authorised persons.

3. Before you start...

Power supply

The scale is powered by batteries. The LC display has a low power consumption – 10,000 weighing operations can be performed with one set of batteries. If  or **bAtE** appear in the display, you should change the batteries. You need 6 1.5 V type AA R6 round-cell batteries.


- Take off the tray and then turn the scale over carefully so that the base is accessible.

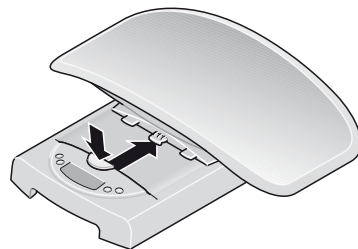- Open the battery compartment lid and insert the batteries supplied according to the polarity shown on the battery compartment.


Check that polarity is correct!

- Now close the battery compartment cover again carefully.

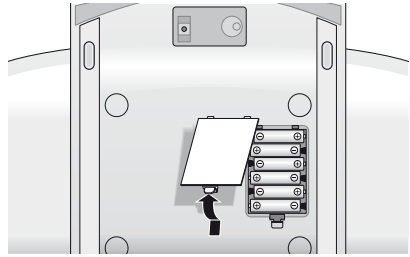
4. Setting up the scale

- Place the scale on a level and stable surface. Then position the tray as shown on the base of the scale and pull forward until the tray is heard to snap into place, and the release button is in the top position. The tray is now securely interlocked with the base.

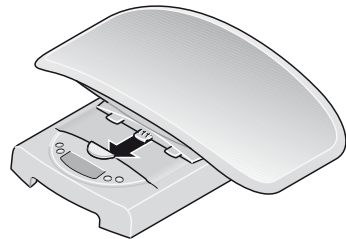- To remove the tray just press the release button hard and slide back the tray.

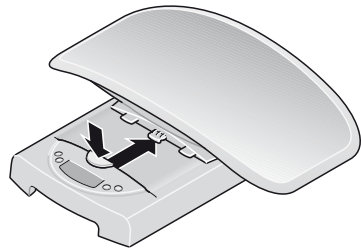5. Operation

Correct weighing


- Press the Start key with no load on the scale.
SECA, **88888** and **0.000** appear consecutively in the display. The scale is then automatically set to zero and ready for use.

- **When weighing infants:** Place the infant on the tray, ensuring that the child lies as still as possible.

- **When weighing children able to stand up:** Remove the tray as described. Check the secure positioning of the base. Place the child on the scale and make sure he/she keeps still.
- Read the weighing result off the digital display.
The scale must not be subjected to weights exceeding the maximum load specified!

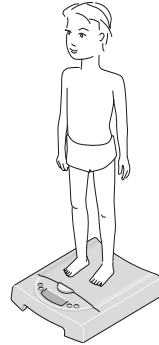- To switch off the scale, press the Start key again.

Automatic switch-off/economy mode

The scale remains ready for weighing for approx. 20 seconds after being switched on. Once weighing has been performed, the value is displayed for approx. 2 minutes until the scale switches off automatically.

Taring weight display (TARE)

The tare function is useful if an additional weight (such as cloths or a pad on the weighing surface) are to be ignored in the weighing result. Proceed as follows:

- Place the additional weight on the scale and press the button **TARE**.

The display of the scale is now reset to "0" and a normal weighing operation can begin.

- You can now perform as many weighing operations as you like with this additional weight.

- TARA is deactivated by pressing the button **TARE** once again, or by switching the scales off.

Important!

The maximum additional weight is 20 kg. The maximum weight which can be displayed is reduced by the weight of the objects already placed on the scale.

Storing the weight value (HOLD)

You can also continue displaying the weight determined once the load has been removed. This makes it possible to put the baby down safely again before noting the weight:

- Press the **HOLD** key whilst the load is on the scale.

The display flashes until a stable weight has been measured. The display is then frozen.

The value remains saved in the display until the hold function is automatically switched off.

- The HOLD function can be switched off by pressing the **HOLD** key again.

0.0 appears in the display if the scale has no load on it, otherwise the current weight. The “HOLD” display vanishes.

Note:

The HOLD and TARE functions can also be used simultaneously.

Use of differential display function

You can use this function for example to measure the intake of milk when breastfeeding or have Baby's increase in weight displayed since being last weighed. This function (which cannot be calibrated) can be used to store an individual weight value and then use that value to form a difference - to detect the amount of milk ingested during breastfeeding, for example, or to see how a particular baby's weight has increased since it was last weighed.

- To store the weight value, press the **STORE** key for approx. 2 seconds after the weighing operation is complete. The **5t0** indicator is faded in. The weight value is now saved and will also remain stored once the scale is switched off.

- With a subsequent weighing operation it is then possible to switch the display over to differential display, i.e. display of current value less saved value, by briefly pressing the **STORE** button.

When the differential display function is activated, **▶◀** appears in the display.

- Briefly pressing the key again switches the function off again.

Note:

If a new weight value is stored, the previous stored value will be overwritten.

If the differential display function is activated with no load on the scale, the stored value is shown in the display as a negative value.

6. Cleaning

Clean the tray and the outside of the scale with a damp cloth as required. Under no circumstances use abrasive or acid cleaners, white spirit, benzene or the like for cleaning.

7. What do I do if...

... no weight is displayed with a load on the scale?

- Check that the batteries are inserted correctly.

... the display shows  ?

- Battery voltage is dropping. Charge the rechargeable batteries in the next few days.

... \overline{bAtt} appears in the display?

- The rechargeable battery is so severely discharged that weighing is no longer possible. Recharge the battery.

... $\overline{5tOP}$ appears in the display?

- Maximum load has been exceeded. Remove the weight on the scale as quickly as possible in order to avoid permanent damage to the scale.

8. Servicing/recalibration

Have the scale recalibrated by authorised persons in accordance with national legal requirements. The year of initial calibration can be found behind the CE symbol on the rating plate above the number of the appointed office, 0109 (Hessen Weights & Measures Office). Recalibration is essential if one or more security marks are damaged.

We recommend having your scale serviced by your local service agent before having it recalibrated. seca Customer Service will be pleased to give you further assistance.

9. Technical data

Dimensions

length:	330 mm
width:	550 mm
height:	165 mm

Weight

approx. 2.5 kg

Temperature range

+10 °C to +40 °C

Maximum load

20 kg

Graduations

20 g

Calibration as per Directive 90/384/EEC

approved for medical use, Class III

Accuracy at first calibration

0 up to 1 kg:	± 10 g
1 kg up to 4 kg:	± 20 g
4 kg up to 20 kg:	± 30 g

Medical product as per Directive 93/42/EEC

Class I

Height of figures

22 mm

Power supply

Battery

Power consumption

approx. 30 mA

10. Disposal

Disposing of the device

Do not dispose of the device in domestic waste. The device must be disposed of properly as electronic waste. Follow the national regulations which apply in your case. For further information, contact our service department at:

service@seca.com

Batteries

Do not throw used batteries away in domestic waste. Dispose of batteries at collection points in the vicinity. When buying new batteries, select those low in harmful substances and containing no mercury (Hg), cadmium (Cd) or lead (Pb).

11. Warranty

A two-year warranty from date of delivery applies to defects attributable to poor materials or workmanship. All moveable parts - batteries, cables, mains units, rechargeable batteries etc. - are excluded. Defects which come under warranty will be made good for the customer at no charge on production of the receipt. No further claims can be entertained. The costs of transport in both directions will be borne by the customer should the equipment be located anywhere other than the customer's premises. In the event of transport damage, claims under warranty can be

honoured only if the complete original packaging was used for any transport and the scale secured and attached in that packaging just as it was when originally packed. All the packaging should therefore be retained.

A claim under warranty will not be honoured if the equipment is opened by persons not expressly authorised by seca to do so.

We would ask our customers abroad to contact their local sales agent in the event of a warranty matter.

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